FORM 5-5

43

Practitioner's Docket No.

PATENT

IN THE UNITED STATES PATENT AND TRADEMARK OFFICE

69.00136

In re application of: Smith et al

Group No.: 1643

Application No.: 09 Filed: 1/31/00

)9 / 494,751

Examiner:

For:

ASSAYS FOR TSH RECEPTOR AUTOANTIBODIES

RECEIVED

APR 19 2000

TECH CENTER 1600/2900

Assistant Commissioner for Patents Washington, D.C. 20231

TRANSMITTAL OF CERTIFIED COPIES

Attached please find the certified copy of the foreign application from which priority is claimed for this case:

Country: United Kingdom

Application Number: 9909661.2

Filing Date: April 28, 1999

Country:

United Kingdom

Application Number: 9812146.0

Filing Date: June 6, 1998

WARNING: "When a document that is required by statute to be certified must be filed, a copy, including a

photocopy or facsimile transmission of the certification is not acceptable 737 GFR. § 1.4(f)

(emphasis added).

SIGNATURE OF PRACTITIONER

Reg. No.: 30,955

KENNETH I. KOHN

(type or print name of practitoner)

Tel. No.: (248) 539-5050

30500 Northwestern Highway, Suite 410

P.O. Address

Farmington Hills, MI 48334

Customer No.:

NOTE: The claim to priority need be in no special form and may be made by the attorney or agent, if the foreign

application is referred to in the oath or declaration, as required by § 1.63.

CERTIFICATE OF MAILING (37 C.F.R. § 1.8a)

I hereby certify that this paper (along with any paper referred to as being attached or enclosed) is being deposited with the United States Postal Service on the date shown below with sufficient postage as first class mail in an envelope addressed to the Assistant Commissioner for Patents, Washington, D.C. 20231.

(type or print name of person mailing paper)

Date: April 11, 2000

Comple Herty

Signature of person mailing paper

(Transmittal of Certified Copies[5-5])

HIS PAGE BLANK (USPTO)

CIBN/103H

and solution of the









The Patent Office
Concept House
Cardiff Road
Newport
South Wales RECEIVED
NP10 8QQ

TECH CENTER 1600/2900

I, the undersigned, being an officer duly authorised in accordance with Section 74(1) and (4) of the Deregulation & Contracting Out Act 1994, to sign and issue certificates on behalf of the Comptroller-General, hereby certify that annexed hereto is a true copy of the documents as originally filed in connection with the patent application identified therein.

In accordance with the Patents (Companies Re-registration) Rules 1982, if a company named in this certificate and any accompanying documents has re-registered under the Companies Act 1980 with the same name as that with which it was registered immediately before re-registration save for the substitution as, or inclusion as, the last part of the name of the words "public limited company" or their equivalents in Welsh, references to the name of the company in this certificate and any accompanying documents shall be treated as references to the name with which it is so re-registered.

In accordance with the rules, the words "public limited company" may be replaced by p.l.c., plc, P.L.C. or PLC.

Re-registration under the Companies Act does not constitute a new legal entity but merely subjects the company to certain additional company law rules.

Signed fursieurs

Dated 24 February 2000

ACE BLANK (USPTO)

atents Act 197.

(Rule 16)

The Patent Office

28APR99 E443988-3 000446..... P01/7700 0.00 - 9909661.1

Request for Grant of a Patent

The Patent Office Cardiff Road NEWPORT NP9 1RH

		NE	WPORT NP9 1RH
1.	Your reference	P79645	
2.	Patent Application number (the Patent Office will fill in this part)	9909661.2	
3.	Full name, address and postcode of the or of each Applicant (underline all surnames)	RSR Limited Avenue Park Pentwyn Cardiff CF2 7HE United Kingdom THE PATENT 2 8 APR 1	r setya.
	Patent Office ADP Number (if you know it)	709425300 RECEIVED B	.
	If the applicant is a corporate body, give the country/state of its incorporation	England and Wales	. 1001
4.	Title of the Invention	Assays for TSH Receptor Autoanti	bodies
5.	Name of your Agent (if you have one)	URQUHART-DYKES & LORD	
	"Address for Service" in the United Kingdom to which all correspondence should be sent (including the postcode)	Alexandra House 1 Alexandra Road SWANSEA SA1 5ED United Kingdom	
	Patents ADP Number (if you know it)	1644005	-
6.	If you are declaring priority from one or more earlier Patent Applications, give the country and the date of filing of the or of each of these	Country Priority application No. (if you know it)	Date of Filing (Day/month/year)
	earlier Applications and (if you know it) the or each	United Kingdom 9812146.0	06.06.1998
7.	If this Application is divided or otherwise derived from an earlier UK Application, give the Number and the Filing Date of the earlier Application	Number of earlier application	Date of Filing (Day/month/year)
8.	Is a Statement of Inventorship and of Right to Grant of a Patent required in support of this request? (Answer 'Yes' if:	Yes	
		•	

a) any Applicant named in part 3 is not an inventor, orb) there is an inventor who is not named as an Applicant, or

c) any named Applicant is a corporate body.)

9.	Enter the number of sheets for any of the following items you are fill with this	
	Form. Do not count copies of the same document	
	Continuation sheet of this Form	0
	Description	15
	Claim(s)	3
	Abstract	1
	Drawing(s)	3 + 3
- 10.	If you are also filing any of the following	
	Priority documents	0
	Translations of priority documents	0
	Statement of Inventorship and Right to Grant a Patent (Patents Form 7/77)	Ö
	Request for Preliminary Examination (Patents Form 9/77)	1
	Request for Substantive Examination (Patents Form 10/77)	0
	Any other documents (please specify)	0
11.		I/We request the grant of a Patent on the basis of this Application
	· .	
		Signature Date
		URQUHART-DYKES & LORD 26 April 1999
12.	Name and daytime telephone number of person to contact in the United Kingdom	URQUHART-DYKES & LORD 26 April 1999 Mr. H.W. Austin - 01792 474327

Warning

After an Application for a Patent has been filed, the Comptroller of the Patent Office will consider whether Publication or communication of the invention should be prohibited or restricted under Section 22 of the Patents Act 1977. You will be informed if it is necessary to prohibit or restrict your invention in this way. Furthermore if you live in the United Kingdom, Section 23 of the Patents Act 1977 stops you from applying for a Patent abroad without first getting written from the Patent Office unless an Application has been filed at least 6 weeks beforehand in the United Kingdom for a Patent for the same invention and either no direction prohibiting Publication or communication has been given, or any such direction has been revoked.

DUPLICATE

Assays for TSH Receptor Autoantibodies

The present invention is concerned with assays (kits and analytical methods) for detecting or monitoring TSH Receptor Autoantibodies.

The hyperthyroidism which is associated with Graves' disease is known to be due to autoantibodies to the thyroidal receptor for thyroid stimulating hormone (TSH). The autoantibodies bind to the receptor and mimic the actions of the natural ligand (TSH), thereby causing the gland to produce high levels of thyroid hormones (as described in Endocrine Reviews 1988, Vol 9, No. 1, pages 106 to 117).

The detection or monitoring of TSH receptor autoantibodies (TRAb) is important in the diagnosis and management of Graves' disease and currently two types of assay are used, namely:

- (a) competitive binding assays which measure the ability of TRAb to inhibit the binding of ¹²⁵I-labelled TSH to preparations of TSH receptor; and
- (b) bioassays which measure the ability of TRAb to stimulate thyroid cells (or other cells transfected with the TSH receptor gene) in culture.

Currently, competitive binding assays (type (a) above) are more widely used, because bioassays of the type mentioned in (b) above are expensive, time-consuming, require highly skilled staff and are unsuitable for wide routine use. In current competitive binding assays, test serum samples $(50\mu l)$ are generally incubated with detergent solubilised porcine TSH receptor $(50 \mu l)$. TRAb present in the test sera bind to the TSH receptor during this incubation. ¹²⁵I-labelled TSH is then added and the incubation continued. During this second incubation, the labelled TSH binds to TSH receptors not already occupied by TRAb. Finally, ¹²⁵I-labelled TSH bound to the receptor is separated from free labelled TSH by addition of polyethylene glycol (PEG), which precipitates the receptor bound TSH but not the free TSH. The radioactivity in the precipitates (separated by centrifugation) is then counted. In the assay, TRAb in test samples inhibits the binding of labelled TSH to the TSH receptor and this results in a lowering of the radioactivity in the precipitates. Assay results can be expressed as an index of inhibition of labelled TSH binding or by use of a set of assay calibrators.

The main limitations of this conventional assay are as follows:

- (a) The assay measures competition between labelled TSH and the TSH receptor and may not detect TRAb which bind well to the receptor but in such a way as not to inhibit TSH binding strongly.
- This results in co-precipitation of all the serum immunoglobulins and the formation of a relatively large pellet. Although the pellet can be counted for radioactivity, it is not a suitable preparation to detect TSH (or other proteins or peptides) labelled with non-radioactive substances such as enzymes or chemiluminescent materials. This is because the serum components in the pellet interfere with such processes as light emission. In addition, the use of PEG precipitation necessitates the use of centrifugation and this is a time-consuming and cumbersome procedure unsuitable for automation.

The present invention concerns procedures and kits for assays for TRAb of the competitive binding type mentioned above. The invention also concerns procedures and kits for TRAb assays of the direct binding type, in which a direct interaction between the receptor and TRAb is used.

According to the present invention, there is provided a method of monitoring autoantibodies to thyroid stimulating hormone (TSH) receptor in a sample of body fluid, comprising, in the following order, the steps of:

- (a) providing (i) porcine TSH receptor or a fragment thereof immobilised to a solid phase, or (ii) TSH receptor which is complexed to a labelled antibody;
- (b) incubating the TSH receptor with a sample of body fluid;
- reacting the incubated sample of body fluid containing the TSH receptor with at least one binding agent which is capable of binding to the TSH receptor in competitive reaction with TSH receptor autoantibodies (TRAb), or when the TSH receptor is (ii), reacting the sample of body fluid, during or after step (b), with at least one binding agent which can bind to TRAb in such a way as not substantially to interfere with binding of the TRAb to the TSH receptor; and
- (d) qualitatively or quantitatively detecting bound TRAb in the reacted incubated sample of body fluid.


According to a preferred feature of the present invention, the sample of body fluid comprises blood, plasma or serum.

The invention comprises the use of antibodies in order to label or immobilise a TSH receptor, the immobilised or labelled TSH receptor being such that it retains its ability to bind TSH and/or TRAb.

The present invention preferably concerns the use of a monoclonal (or polyclonal) antibody to the TSH receptor, which is bound strongly to the receptor at a site distinct from the part of the receptor which binds TSH and TRAb. The antibody can bind to the receptor strongly at the same time as TSH or TRAb and can be used to alleviate many of the limitations of the current TRAb assay method.

In a first embodiment of the invention, the antibody can be immobilised on a solid phase (such as a plastic tube or plastic plate, or magnetic or non-magnetic particles) using standard procedures. This solid phase can then be used instead of PEG to separate labelled TSH bound to the TSH receptor from free labelled TSH. The TSH (or similar ligand) can be labelled with isotopic or non-isotopic labels.

There are several applications of this first embodiment, of which the following are schematic examples:

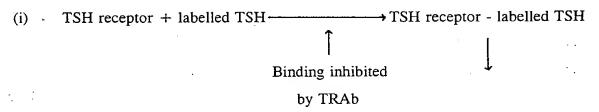

separate from unbound labelled TSH by addition of immobilised receptor antibody (□-Ab)

□-Ab-TSH receptor-labelled TSH

(ii)	□-Ab-TSH receptor + labelle	ed TSH ———→ □-Ab-TSH receptor-labelled TSF
		↑
	Form this complex first	binding inhibited
		by TRAb

(iii) In (i) and (ii) above, the receptor Ab can be bound **indirectly** to the solid phase. For example, the Ab can be biotinylated and reacted with the solid phase which contains avidin or streptavidin:

(a) Bi-Ab-TSH receptor + labelled TSH → Bi-Ab-TSH receptor-labelled TSH complex

binding inhibited by TRAb separate from unbound labelled TSH by addition of immobilised avidin or streptavidin (□-Avidin)

□-Avidin-Bi-Ab-TSH receptor-labelled TSH

In a second embodiment of the invention, the antibody can be labelled directly with an isotopic label such as ¹²⁵I, or with a non-isotopic label such as an enzyme, dye, or chemiluminescent compound. Alternatively, the antibody may be labelled indirectly using, for example, the avidin-biotin system. The labelled antibody can then be used to label the TSH receptor itself and this complex of receptor and antibody can then be used to detect and/or monitor TRAb (either qualitatively or quantitatively).

Examples of this approach are described schematically as follows:

(i) Labelled Ab - TSH receptor + TRAb → labelled Ab-TSH receptor

Form this complex first

TRAb

Add anti-human IgG reagent such as solid phase protein A (□-protein A)

Labelled Ab-TSH receptor

□ - protein A - TRAb

In this system, increasing amounts of TRAb in a test sample will result in increasing amounts of labelled Ab in the final complex.

(ii) This example depends on the coupling of TSH (or similar ligand) directly or indirectly to a solid phase (□-TSH) receptor labelled with the antibody (which is itself labelled with an isotopic or non-isotopic label), which can bind to this immobilised TSH (or similar ligand). This binding will be inhibited by TRAb, as illustrated below:-

The schemes described above with reference to the first and second embodiments of the invention illustrate certain advantages of the present invention. The immobilised antibody may be used to label or immobilise the TSH receptor in such a way that the receptor retains its ability to bind TSH and/or TRAb. Alternatively, the immobilised antibody may be used for the monitoring of TRAb (or other ligands which interact with the TSH receptor) in patient sera, or for purification of the TSH receptor.

In particular, the ability to label the TSH receptor using the antibody permits monitoring of the direct interaction of TRAb with the TSH receptor as in the second embodiment of the invention described above. Furthermore, the present invention allows the immobilisation of the receptor before, during or after interaction with TSH and TRAb. This ability to immobilise the receptor can be used to create new TRAb assays which do not depend on PEG and/or isotopic labels. Such assays may be suitable for automation and immunochromatographic systems.

The TSH receptor (TSHR) is present in very low numbers on the surface of thyrocytes (about 10³ receptors per cell) which has made the receptor very difficult to purify from native sources (as described in "Baillière's Clinical Endocrinology and Metabolism", 1997, Vol II, pages 451 to 474 - Sanders).

In contrast, recombinant TSHR can be expressed in mammalian cells (for example, in Chinese hamster ovary (CHO) cells) at much higher levels of about 10⁵-10⁶ receptors per cell (Sanders). In addition, recombinant TSHR preparations produced in non-thyroid cells are not contaminated with other thyroid autoantigens such as thyroglobulin or thyroid peroxidase (Springer Seminars in Immunopathology, 1993, pages 309-318 - Furmaniak).

Recombinant TSHR preparations produced in mammalian cells are the only ones which show TSH and TRAb binding characteristics similar to that of the native receptors. Such binding characteristics are not, for example, produced in yeast, insect cells or bacteria. This is because of a very complex relationship between the TSHR's structure and the TSH/TRAb binding sites. The TSHR's post-translational processing and folding of the "mature" protein is best achieved in the mammalian cell environment (see Sanders).

Purification of large amounts of recombinant TSHR from mammalian cells with its TSH and TRAb binding activities intact has not been reported. One of the major problems is loss of TSH/TRAb binding activity following binding with mouse monoclonal antibodies interacting with the extracellular part of the TSHR (Sanders). Also, to date, development of new and convenient strategies for routine measurement of TRAb has not been successful.

Preferred features of the present invention are illustrated by the following nonlimiting detailed worked examples.

EXAMPLES

1. Cloning of porcine TSHR cDNA

RNA was extracted from porcine thyroid tissue using the guanidinium thiocyanate-phenol chloroform method (Chomczynski and Sacchi, Anal. Biochem. Vol 162. 1987, pages 156 to 159). mRNA was purified from total RNA using a Dynal bead mRNA purification kit (Dynal, Wirral L62 SAZ UK). This mRNA was used to make a cDNA library using ZapExpress cDNA Gigapack Cloning Kit III (Strategene Ltd., Cambridge CB4 4DF UK). Four degenerate oligonucleotides were made to the known TSHR sequences (mouse, rat, human, dog and bovine) and two fragments of porcine TSHR amplified using PCR. These were sequenced to verify their homology with TSHR cDNA and used to screen the cDNA library for full-length porcine TSHR clones. Three full-length clones were obtained and fully sequenced.

2. Expression of porcine recombinant TSHR protein in CHO cells

An ATG start codon in the 5' untranslated region (5'UTR) of the full-length pTSHR cDNA was removed by PCR and the cDNA cloned into pcDNA 3.1 (+) (Invitrogen BV, 9351 NV Leek, The Netherlands). DNA encoding the full-length TSHR was transfected into CHO cells (CHO-KI from ECACC, Porton Down SP4 OJG UK) by electroporation. Clones expressing TSHR were detected using ¹²⁵I-TSH binding directly to cells growing on 24-well plates. The clones showing highest TSH binding were expanded and recloned twice by limiting dilution. One stable cell line expressing about 4 x 10⁵ TSHR per cell was chosen for expansion and production of recombinant TSHR.

3. Preparation of detergent solubilised recombinant porcine TSHR

CHO cells expressing TSHR were grown to confluence, detached from roller bottles and the cell pellet washed with ice-cold 50 mM NaCl, 10 mM Tris-HCl, pH 7.5 containing 1 mM phenylmethylsulfonylfluoride (PMSF), then homogenised in the same buffer. The cell membranes after centrifugation at 12000g for 30 min at 4°C were solubilised in the same buffer (4 ml of buffer for approximately 4 x 10⁸ cells) as used for homogenisation except for addition of 1% Triton X-100. The solubilised receptor preparations were centrifuged at 90,000g for two hours at 4°C and the supernatant stored at -70°C.

4. Expression of the C-terminal end of the porcine TSHR protein

The expression in Escherichia coli as a fusion protein with glutathione S-transferase (GST) was carried out using standard protocols (as described in Journal of Molecular Endocrinology (1998) Vol 20, pages 233-244 - Oda). The 3' end of cDNA (1809 to 2295 bp) coding for the last 160 amino acids was cloned in frame with the GST fusion protein in pGEX2T vector (Pharmacia Biotech, St. Albans ALl 3AW UK). An overnight culture of E. coli (strain UT580) transformed with pGEX-2T/TSHR plasmids was diluted 1/5 into 2 x YTG medium (16 g Tryptone, 10g yeast extract, 5 g NaCl, 20 g glucose per litre, pH 7.0), incubated for 3 hours at 30°C. Thereafter, isopropyl-3-D-thiogalactopyranoside (IPTG) was added to a final concentration of 1 mM in order to induce protein expression, followed by incubation for a further three hours. The bacterial pellets were resuspended in PBS (8 g NaCl, 0.2 g KCl, 1.44 g Na₂HPO₄ and 0.24 g KH₂PO₄ per litre, pH 7.4) containing 1% Triton X-100 and sonicated three times for one minute on ice. The inclusion bodies were pelleted, washed in 4M urea, solubilised in 8M urea and separated on 9% polyacrylamide gels (SDS-polyacrylamide electrophoresis, SDS-PAGE) under reducing conditions. The TSHR/GST fusion proteins (mol. wt. 44 kDa) were electroeluted from polyacrylamide gel slices in 0.1M NaHCO3 and 0.1% SDS pH 7.8, dialysed against 50 mM Tris-HCl pH 8.0 and stored in aliquots at -70°C.

5. Preparation and purification of monoclonal antibodies to the porcine TSHR

Electroeluted TSHR/GST protein was used to immunise BALB C mice (50 μg per mouse per injection) until the titre of antibody to the TSHR was high. The TSHR antibody level in mouse sera was tested using immunoprecipitation assay based on ³⁵S-labelled TSHR produced in an in vitro transcription/translation system (see the method described in Journal of Clinical Endocrinology and Metabolism, Vol 82 (1997) No. 4, pages 1288-1292 - Prentice).

Mouse spleen cells were fused with mouse myeloma cell line (X63-Ag8.653 from ECACC) and cloned to produce stable hybridomas secreting the TSHR antibody using standard techniques (Oda). The antibody 4E31 was found to precipitate ³⁵S-TSHR in the immunoprecipitation assay and to react well with the TSHR in western blotting analysis.

4E31 was purified from hybridoma culture supernatants by chromatography on a Prosep-A (Bioprocessing, Consett DH8 6TJ UK) column. Fab₂ fragments were obtained following digestion with pepsin (pepsin concentration of 1 mg/ml in 70 mM sodium acetate, 50 mM NaCl pH 4.0, IgG to enzyme ratio of 5:1) and chromatography on Prosep A to adsorb intact IgG. 4E31 Fab₂ preparations of about 1 mg/ml were stored at -70°C in aliquots.

6: TSH receptor binding characteristics of the 4E31 antibody

Solubilised recombinant porcine TSH receptor was incubated with 125 I-labelled TSH to form a 125 I-labelled TSH-TSH receptor complex. The 4E31 IgG (or control IgG) was immobilised by linking it to magnetic latex beads and these beads (100 μ l; 1 μ g of 4E31) were incubated with the preformed 125 I-labelled TSH-TSH receptor complex (100 μ l). After 1 hr at 37°C, the beads were separated on a magnet, washed and counted for 125 I. A sample of the 125 I-labelled TSH-TSH receptor complex was precipitated with PEG to determine the amount of free labelled TSH present. Table 1 shows the results obtained:

Table 1

Sample	cpm Bound to Beads
4E31 IgG beads plus:-	
(a) labelled TSH-TSH receptor complex	16,374
(b) labelled TSH only	2,118
Control IgG beads plus labelled TSH-TSH receptor complex	2,433
Labelled TSH-TSH receptor complex precipitated by PEG	16,499
Total cpm in 100 μ l of labelled complex	37,073

(Control IgG used was a mouse monoclonal antibody to glutamic acid decarboxylase)

These studies indicated clearly that 4E31 could bind to the TSH receptor at the same time as TSH. Further, 4E31 coupled to a solid support could be used to separate TSH bound to the TSH receptor from free TSH with results similar to those obtained with PEG.

Table 2 shows that in the presence of individual healthy human donor sera, about 55 % of the labelled TSH added bound to tubes coated with 4E31 and recombinant porcine TSH receptor (according to the first embodiment of the invention, scheme ii above). Similar binding was observed in the case of sera from patients with Hashimoto's thyroiditis and patients with systemic lupus erythematosus (SLE). However with sera from the 5 patients with Graves' disease, there was markedly less labelled TSH binding (about 35 %). All 5 sera contained readily detectable amounts of TSH receptor autoantibody as judged by inhibition of TSH binding to native porcine TSH receptor and separation of receptor bound and free TSH with PEG.

Table 2 Binding of labelled TSH to plastic tubes coated with 4E31 (Fab)₂ followed by recombinant porcine TSH receptor - effect of different sera.

		inhibition of TSH binding (%)		
Sample	cpm bound	coated tube method (recombinant receptor)	PEG method (native receptor)	
Individual normal				
sera:-				
1	15,171		7.74	
2	15,209	1		
3	15,480		A. Carrier	
4	14,768			
5	15,496			
Individual Graves'				
sera:-			· A.	
1	10,132	34	27	
2	9,575	37	31	
3	12,495	18	21	
4	11,478	25	11	
5	9,163	40	. 48	
Individual				
Hashimoto sera:-			.\$	
1	14,659			
2	15,215			
3 .	14,603			
4	15,026			
5	15,370			
Individual SLE				
sera:-				
1	15,711			
2	15,261			
3	15,091			
4	15,386			
5	15,403			
Total cpm	27,420			

Table 3 shows similar results with tubes coated with native (i.e. non-recombinant) porcine TSH receptor via 4E31. Again, the results obtained with the coated tube assay are similar to those obtained with the PEG separation method.

Table 3 Binding of labelled TSH to plastic tubes coated with 4E31 (Fab)₂ followed by native detergent solubilised porcine TSH receptor.

	inhibition of TSH binding (%) with:-		
Sample	native porcine TSH receptor bound to plastic tubes (via 4E31 (Fab) ₂)	native porcine TSH receptor PEG precipitation assay	
Individual normal sera:-			
· a	1.5	-0.62	
b	2.2	3.2	
С	-5.0	3.9	
d	1.8	, 3.2	
Graves sera:-			
6	74	71	
7	75	72	
8	81	76	
9	70	73	
10	77	88	
11	88	83	
12	88	72	
13	80	76	
14	57	70	
15	90	86	
16	55	72	
17	81	68	

The inhibiting effect on ¹²⁵I-TSH binding shown in tables 2 and 3 was dependent on the concentration of TRAb in the serum. As shown in Table 4, increasing amounts of TRAb standard preparation (thyroid stimulating antibody 1st International Standard 90/672) had a dose dependent effect on TSH binding inhibition.

Table 4

TRAb standard dilution (mu/ml)	native porcine TSH receptor (PEG separation method) - inhibition of TSH binding (%)	4E31/recombinant porcine TSH receptor coated tubes - inhibiton of TSH binding (%)
0	0	0 .
1.25	5.1	5.1
2.5	9.9	12.0
5	16.2	17.0
10	28.4	31.7
20 .	44.5	51.6
40	69.1	74.8

Direct precipitation of ¹²⁵I-4E31 Fab₂-TSH complexes

Another example of the invention (according to the second embodiment of the invention, scheme (i) above) is an assay in which 4E31 (Fab)₂ is labelled with ¹²⁵I and reacted with the TSH receptor. This complex is then incubated with TSH receptor autoantibodies in patients' sera and the resulting termolecular complex precipitated by addition of solid phase protein A. An example is shown in Table 5.

Table 5 Reaction between recombinant porcine TSH receptor labelled with $^{125}\text{I-}4E31$ (Fab)₂ and TSH receptor autoantibodies.

Sample	cpm bound
Labelled receptor plus:-	
Individual normal serum:-	
e f	1,994 2,041
90/672 reference material (mU/ml): 1.25 2.5 5 10 20 40	2,715 2,958 3,358 4,202 5,868 8,949
Graves' serum:- 18 19	12,326 13,225
Total cpm	48,428

Graves' sera 18 and 19 gave TSH binding inhibition values of 31% and 73% respectively in the PEG precipitation method using native porcine TSH receptor.

4E31 could be labelled with biotin, bound to streptavidin coated tubes and then reacted with porcine TSH receptor - see the first embodiment of the invention, Scheme (iii) b. Receptor immobilised in this way readily bound ¹²⁵I-labelled TSH and this binding was inhibited by TSH receptor autoantibodies in patient sera. An example is shown in Table 6.

Table 6 Labelled TSH binding to TSH receptor - bound to streptavidin coated tubes via biotinylated 4E31 (Fab)₂ and effect of TRAb.

Serum Sample	125I-labelled TSH bound (cpm)
Individual healthy normal sera:-	
g	. 12,225
h	12,164
TRAb positive Graves' sera:-	
20	2,971
21	7,400
Total cpm	27,339

The present invention will now be illustrated with reference to the following figures which are by way of example only.

Figure 1 illustrates the correlation between results using a method of monitoring autoantibodies to thyroid stimulating hormone receptor using TRAb coated tubes and a method using PEG precipitation assays, both methods using ¹²⁵I-labelled bovine TSH. Furthermore, results in different patient groups with the method using TRAb coated tubes are shown.

Figure 2 illustrates results of an ELISA method using TSH receptor coated plates; 2a shows the effect of TRAb standard 90/672 on porcine TSH binding to receptor-coated plates; 2b shows a comparison of an ELISA method and a PEG precipitation assay; and results in different patient groups with a TRAb ELISA method using bovine TSH are shown in 2c.

Figure 3 illustrates the correlation between a conventionally used TRAb assay method and a direct precipitation assay method according to the invention. Furthermore, results in different patient groups with a direct TRAb assay are shown.

Claims:

- 1. Method of monitoring autoantibodies to thyroid stimulating hormone (TSH) receptor in a sample of body fluid, comprising, in the following order, the steps of:
 - (a) providing
 - (i) porcine TSH receptor or a fragment thereof immobilised to a solid phase; or
 - (ii) TSH receptor which is complexed to a labelled antibody;
 - (b) incubating said TSH receptor with a sample of body fluid;
 - (c) reacting said incubated sample of body fluid containing said TSH receptor with at least one binding agent which is capable of binding to said TSH receptor in competitive reaction with TSH receptor autoantibodies (TRAb), or when said TSH receptor is (ii), reacting said incubated sample of body fluid, during or after step (b), with at least one binding agent which can bind to TRAb in such a way as not substantially to interfere with binding of said TRAb to said TSH receptor; and
 - (d) qualitatively or quantitatively detecting bound TRAb in said reacted incubated sample of body fluid.
- 2. Method according to claim 1, wherein said sample of body fluid comprises blood, plasma or serum.
- 3. Method according to claim 1 or 2, wherein said solid phase comprises a plastics material, a magnetic material or a non-magnetic material.
- 4. Method according to any of claims 1 to 3, wherein said labelled antibody comprises an antibody to TSH receptor.
- 5. Method according to any of claims 1 to 4, wherein said agent capable of binding to said TSH receptor comprises any of TRAb, mouse monoclonal antibodies, human monoclonal antibodies, peptides and recombinant antibodies.
- 6. Method according to any of claims 1 to 5, wherein said agent capable of binding to said TSH receptor is immobilised to a solid phase.

- 7. Method according to any of claims 1 to 6, wherein said agent capable of binding to said TSH receptor is labelled isotopically.
- 8. Method according to claim 7, wherein said isotopic label comprises ¹²⁵I.
- 9. Method according to any of claims 1 to 6, wherein said agent is labelled non-isotopically by means of an enzyme, dye, or fluorescent or chemiluminescent material.
- 10. Method according to any of claims 1 to 9, wherein said porcine TSH receptor is immobilised, either directly or indirectly, to said solid phase, via an antibody to the TSH receptor.
- Method according to any of claims 1 to 10, wherein said labelled antibody is indirectly labelled with an organic compound.
- 12. Method according to claim 11, wherein said organic compound comprises biotin.
- 13. Method according to claim 11 or 12, wherein said organic compound is complexed to a protein with an affinity for said compound.
- 14. Method according to claim 13, wherein said protein comprises avidin or streptavidin.
- 15. Method according to any of claims 1 to 14, wherein said porcine TSH receptor is derived from porcine thyroid tissue.
- 16. Method according to any of claims 1 to 14, wherein said porcine TSH receptor comprises recombinant material with at least one epitope for TRAb.
- 17. Method according to any of claims 1 to 14, wherein said TSH receptor comprises recombinant or synthetic material with at least one epitope for TRAb.
- 18. Method according to any of claims 1 to 17, wherein step (d) comprises detecting of labelled or unlabelled TSH bound to TSH receptor and unbound labelled or unlabelled TSH.

- 19. Method according to claim 18, wherein said labelled or unlabelled TSH comprises bovine or porcine TSH.
- 20. Method according to any of claims 1 to 17, wherein step (d) comprises detecting of labelled or unlabelled TSH agonist bound to TSH receptor and unbound labelled or unlabelled TSH agonist.
- 21. Method according to any of claims 1 to 17, wherein step (d) comprises detecting of labelled or unlabelled TSH antagonist bound to TSH receptor and unbound labelled or unlabelled TSH antagonist.
- 22. Method according to claim 20 or 21, wherein said TSH agonist or TSH antagonist is a monoclonal antibody reactive with said TSH receptor, or a fragment of said monoclonal antibody.
- 23. Method according to claim 22, in which said monoclonal antibody is a human monoclonal antibody.
- 24. Method according to claim 22, wherein said monoclonal antibody is a recombinant antibody or recombinant antibody fragment.
- 25. Method according to claim 20 or 21, wherein said TSH agonist or TSH antagonist is a peptide.
- 26. Method according to claim 25, wherein said peptide is derived from TSH.

ABSTRACT

Assays for TSH Receptor Autoantibodies


A method of monitoring autoantibodies to thyroid stimulating hormone (TSH) receptor in a sample of body fluid, comprising the steps of:

- (a) incubating TSH receptor with a sample of body fluid;
- (b) reacting the incubated sample of body fluid with at least one binding agent which is capable of binding to the TSH receptor in competitive reaction with TSH receptor autoantibodies (TRAb), or in a case where TSH receptor is complexed to labelled antibody, reacting the sample of body fluid with at least one binding agent which can bind to TRAb in such a way as not substantially to interfere with binding of the TRAb to the TSH receptor; and
- (c) detecting bound TRAb in the reacted incubated sample of body fluid.


THIS PAGE BLANK (USPTO)

Figure 1

COMPARISON OF TRAB COATED TUBES (CT) AND PEG PRECIPITATION ASSAYS - BOTH USING 125-I-LABELLED TSH

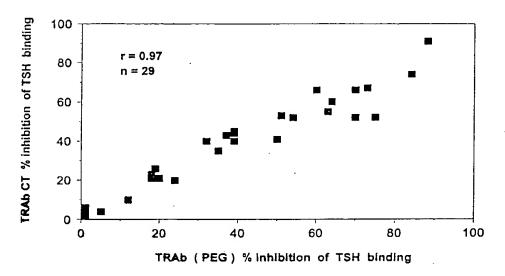


TRAb COATED TUBE ASSAY - RESULTS IN DIFFERENT PATIENT GROUPS

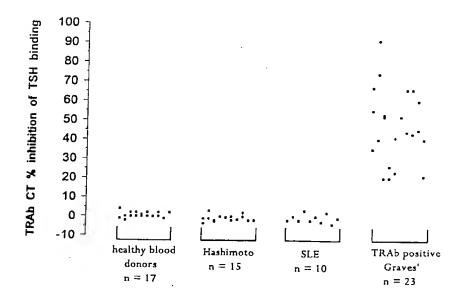


THIS PAGE BLANK (USPTO)

Figure 2


TRAB ELISA USING TSH RECEPTOR COATED PLATES INHIBITION OF PORCINE TSH-PEROXIDASE BINDING BY TRAB STANDARD 90/672

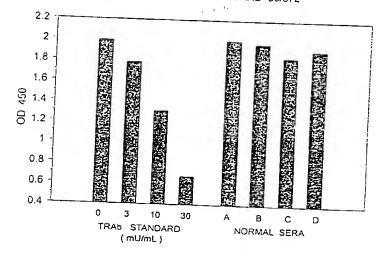

COMPARISON OF ELISA AND PEG PRECIPITATION ASSAYS FOR TRAB USING 21 GRAVES SERA

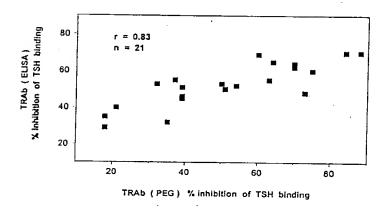



TRAD ELISA RESULTS IN DIFFERENT PATIENT GROUPS

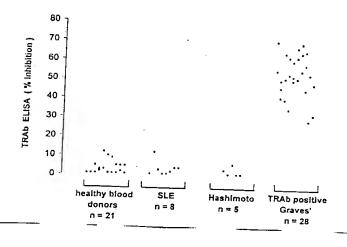


THIS PAGE BLANK (USPTO)

Figure 3

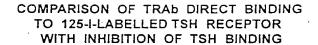

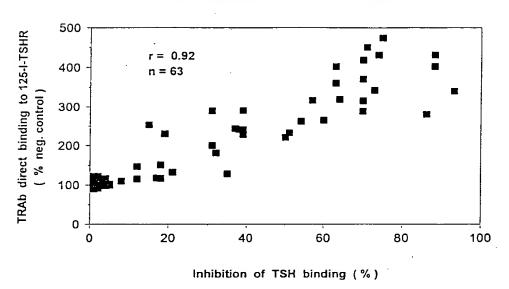

DIRECT TRAB ASSAY - RESULTS IN DIFFERENT PATIENT GROUPS

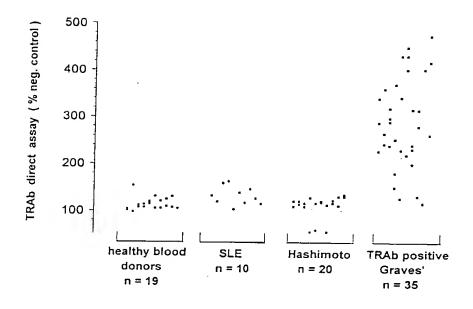

THIS PAGE BLANK (USPTO)